



WIRELESS CONNECTIVITY

ReSound Micro Mic Quick guide

– GETTING STARTED



resound.com

B STATUS LIGHT


- Will blink green every other second when the Micro Mic is turned on and charged.
- Will blink orange every second when there is less than 30 minutes of battery left.

C VOLUME CONTROL

Use the volume buttons on the side of the Micro Mic to adjust the volume. Press '+' to increase the volume and '-' to decrease the volume. The Mute/Volume light **E** will blink each time the volume is adjusted.

D MUTE

To mute the Micro Mic press the mute button. Now the mute/volume status light **E** will blink orange every second. To unmute press the mute button again, the mute/volume light will stop blinking.

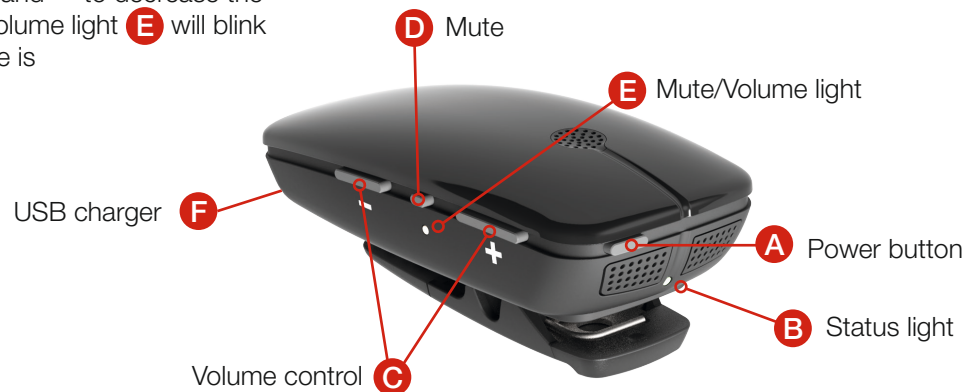

GETTING STARTED

(Charge the Micro Mic before use)

1. Turn on the Micro Mic: Push and **hold** the power button **A** on the top of the Micro Mic until status light **B** blinks green.
2. **Connect your hearing aids** to the Micro Mic: Push and **hold** the program button/toggle button **on your hearing aids** until you hear a melody. Check the connection by blowing or tapping on the Micro Mic. This sound should be audible in your hearing aids. If this is not the case, see the troubleshooting guide on the reverse page.
3. To exit the Micro Mic streaming program: Push the program button/toggle button on your hearing aids. Your hearing aids will return to the most recent hearing aid program indicated by one or more beeps.
4. Turn off the Micro Mic: Push and **hold** the power button until the status light turns red and then off.



F CHARGING

- Insert the charger in the micro USB plug in the bottom of the Micro Mic.
- Status light indications when charging:



Power/Charging	Charging	Fully charged
Powered Off		
Powered On		

TO PAIR YOUR HEARING AIDS AND MICRO MIC

(Your dispenser will usually do this for you)

1. Push the power button on the Micro Mic and hold until the status light turns green. The Micro Mic is now turned on.
2. Find the pairing button placed below the pairing icon (📶) on the back side of the Micro Mic. If you cannot find this icon, try rotating the clip.
3. Push the pairing button (use a pen or similar). The status light will now blink orange.
4. Open and close the battery doors on the hearing aids. This must be done within 20 seconds of pushing the pairing button.


A successful pairing will be indicated by a melody playing in the hearing aids and the status light will stop blinking and be solid orange for 2 seconds before resuming blinking.



Note: Remember that the hearing aids need to be switched into the Micro Mic streaming program to listen to the Micro Mic

The Micro Mic sends speech directly to your hearing aids. It can be clipped to the desired talker's clothes, hung on a lanyard or be placed on the table.

Speech understanding is easier when using the Micro Mic because the microphone is closer to the speaker's mouth than your hearing aids are.

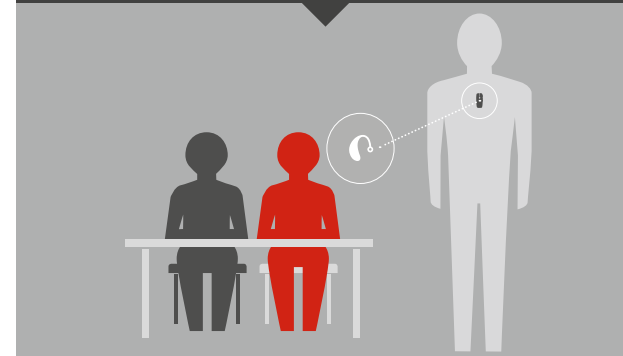

TROUBLESHOOTING GUIDE – NO CONNECTION

1. **Check the status light.** It should blink green every 2 seconds. If it is not blinking, try pushing and holding the power button until you see a green light. If the status light turns red and off the Micro Mic needs to be charged. If the power status light blinks orange, this indicates that the Micro Mic needs to be charged.
2. **Check the Mute/Volume light,** if it is blinking orange this indicates that the Micro Mic is muted. No sound is being sent to the hearing aids. Push the mute button to unmute.
3. **Check hearing aid battery status:** Change the battery in your hearing aids. Hearing aids sometimes need new batteries to receive a streamed signal, even if they are fine for normal hearing aid use.
4. **If Multiple Streaming Devices are connected:** If you have a TV streamer or more than one Micro Mic
 - Try pushing the hearing aid program button/toggle button and **holding** until you hear a melody
 - Repeat this step (two long presses) (Micro Mic can be placed in 'channel' 1, 2 or 3.)
5. **Re-pair the Micro Mic to the hearing aids:** Directions on how to pair your hearing aids and the Micro Mic can be found above.

ReSound North America

8001 Bloomington Freeway
Bloomington, MN 55420
1-888-735-4327
resound.com

ReSound Canada

303 Supertest Road
Toronto, Ontario M3J 2M4
1-888-737-6863
resound.com

ReSound Government Services

8001 Bloomington Freeway
Bloomington, MN 55420
1-800-392-9932
resound.com/veterans

For further instructions, please refer to the full user guide provided with your Micro Mic